Portable Battery Pack

User Guide / Guide de l'utilisateur / Guida dell'utente /
Gebrauchsanweisung / Guía del usuario / Guia do Utilizador /
Guia do Usuário / Οδηγίες χρήσης / Betjeningsvejledning /
Gebruikershandleiding / Användarhandbok / Käyttöohje /
Bruksanvisning / Instrukcja obsługi / Uživatelská příručka /
Használati utasítás / Používateľská príručka / Navodila za uporabnike /
Lietotāja rokasgrāmata / 使用者指南 / Руководство пользователя /
Kullanım Kılavuzu / 사용자 설명서 / ユーザーガイド

Portable Battery Pack User Guide

Introduction

The Portable Battery Pack is an accessory for use with the DC Cord or DC Power Adapter to operate many Respironics therapy devices from a 14.4 AH battery.

Package Contents

The items to the right are included with your Portable Battery Pack. If any of the items shown are missing, contact your home care provider.

Battery Charger

Carrying Case with
Battery

USER
GUIDE

Power Cord

NOTE: If you ordered the Portable Battery Pack with the DC Power Adapter or DC Cord, additional items will be included in your package. If you ordered a replacement carrying case with the battery, you will not receive a charger.

1068511
DSF 08/28/2009

English

Warnings and Cautions

Warnings

- Do not store the battery pack in a hot environment. Storage in a hot environment (e.g., automobile in hot weather, under direct sunlight, in front of a stove, or near a fire) may corrode the battery, ignite, or create an electrical shock hazard.
- If the battery leaks and the gel spills on the skin or clothing, immediately wash it off with plenty of water. If the gel comes in contact with your eyes, immediately flush your eyes with plenty of clean water and seek medical attention.
- The batteries should be used in a ventilated area.
- Do not underestimate the weight of the battery.

Cautions

- Avoid dropping the battery pack. Damage and deterioration of the battery may occur if the battery is dropped.
- Replace the battery (PN 1028871) if it fails to charge or if you notice a decrease in capacity.
- Inspect the power cord for signs of damage. Discontinue use and replace the charger (enclosure or cord) if it is damaged.

Refer to the instructions included with your DC Cord or DC Power Adapter for additional warnings and cautions.

Charging the Battery

IMPORTANT! The battery must be in an upright position during use and while charging.

IMPORTANT! Remove the charger from the carrying case before charging.

1. Connect the auto adapter end of the charger into the battery pack socket, located in the front pouch of the carrying case. See Figure 1.
2. Plug the power cord into the charger.
3. Monitor the LED on the charger while plugging the power cord into the wall outlet and ensure that it turns RED (the charger LED may turn RED for only a short time before turning ORANGE). If the charger LED is RED for any length of time, the charger is correctly connected to the battery. If the charger LED does not turn RED for any length of time, leave the power cord disconnected from the wall outlet for 3 minutes and repeat this step. If the charger LED does not turn RED, ensure that the auto adapter is fully seated in the battery pack socket and that the battery and charger fuses are not blown (see Replacing the Battery Fuse and Replacing the Charger Fuse sections). If the charger LED still does not turn RED, contact your home care provider. During normal operation, the battery charger's LED will appear RED while it is charging. When the battery is 80-95% charged, the LED will appear ORANGE. When the charger is fully charged, the LED will turn GREEN.
4. After the LED turns GREEN, the charger can be disconnected from the battery.

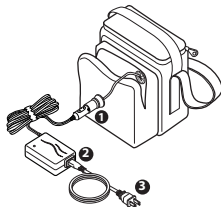
Figure 1: Charging the Battery

Replacing the Battery Fuse

If the battery fails to charge or operate your device, check the fuse and replace it if it has blown, using the instructions below. If the problem continues, contact your home care provider.

1. Pull the battery cable through the opening in the front pouch of the carrying case until the fuse holder is accessible.
2. Remove the fuse holder. See Figure 2.
3. Replace the fuse with a standard mini 15 A blue fuse.
NOTE: To ensure proper operation of your battery, use a 15 A fuse only.
4. Place the fuse holder back onto the fuse.

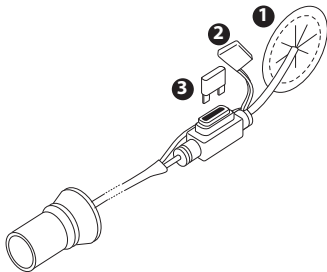
Figure 2: Replacing the Battery Fuse

Replacing the Charger Fuse

If the battery fails to charge, check the fuse, located in the auto adapter. If the fuse has blown, replace it using the instructions below. If the problem continues, contact your home care provider.

1. Open the fuse holder by turning the retaining nut counter-clockwise.
2. Remove the tip of the fuse holder.
3. Replace the fuse with a 3 A fuse, being careful not to dislodge the spring inside the fuse holder.
4. Replace the tip of the fuse holder and retaining nut. Turn the retaining nut clockwise until the connection is tight.

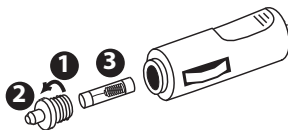
Figure 3: Replacing the Charger Fuse

IMPORTANT! The battery must be in an upright position during use and while charging. Actual power consumption will vary according to prescription pressures and breathing patterns. The lower the prescription, the less power required.

NOTE: The battery should be fully charged before use to ensure that the full battery capacity is available.

NOTE: If you are using a heated humidifier with your therapy device, the humidifier will not heat up while connected to DC power. In this mode, the humidifier acts as a passover, or non-heated humidifier.

NOTE: The battery adapter cable accessory, used with the DC Power Adapter, cannot be used with the Portable Battery Pack.

Connecting the DC Power Adapter

1. Turn off the therapy device.
2. Connect the Auto Adapter to the input cable end of the DC Power Adapter. See Figure 4.
3. Plug the Auto Adapter into the battery pack socket.
4. Plug the DC Power Adapter's output cable into one of the power receptacles on the rear of the therapy device.
5. Turn therapy device on, and verify proper operation.

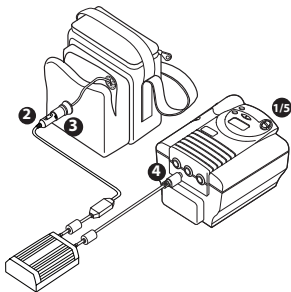
Figure 4: Connecting the DC Power Adapter

Connecting the DC Cord

1. Turn off the therapy device.
2. Plug the adapter end of the DC Cord into the battery pack socket. See Figure 5.
3. Plug the connector end of the DC Cord into the therapy device. The LED on the DC Cord will light when the connections are complete.
4. Turn the therapy device on, and verify proper operation.

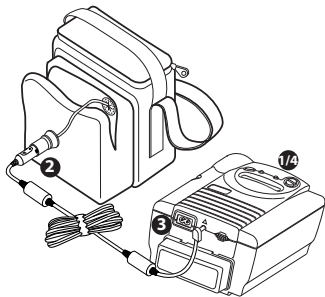
Figure 5: Connecting the DC Cord

Specifications

Battery Pack

Power:	14.4 AH	Actual power consumption will vary according to prescription pressures and breathing patterns. The lower the prescription, the less power required.	
Weight:	6.2 kg		
Fuse:	Mini 15 A (Littlefuse PN 297015, Respironics PN 1044137 or equivalent)		
Conditions:		Operating	Storage
	Temperature:	0°C to 40°C (charge) -20°C to 50°C (discharge)	-15°C to 40°C
	Relative Humidity:	25-85%	25-85%

Charger

Input Rating:	Nominal 90-264 VAC/47-63 Hz
Nominal Voltage Rating:	12 V
Charging Voltage:	14.7 V
Charging Current:	2.3 A
Auto Adapter Fuse:	3 A
Maximum Output Power:	35 W
Temperature Range:	-10°C to 40°C
Safety:	Conforms to EN 60601-1

Disposal

Symbol	Description
	Compliant with the Waste Electrical and Electronic Equipment and Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Recycling Directives (WEEE/RoHS).
	The battery is fully recyclable and should be accepted at any location that accepts automotive starter batteries.

Reordering

Reorder Number	Description
1028869	Portable Battery Pack (includes carrying case with battery, battery charger, and power cord)
1028870	Battery Charger
1028871	Battery with Carrying Case

For additional information, contact Respiration Customer Service at 1-724-387-4000 or Respiration International Customer Service at +49 8152 93060.

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Respironics Inc.